


# Disposable Hemocytometer

# C-CHiP

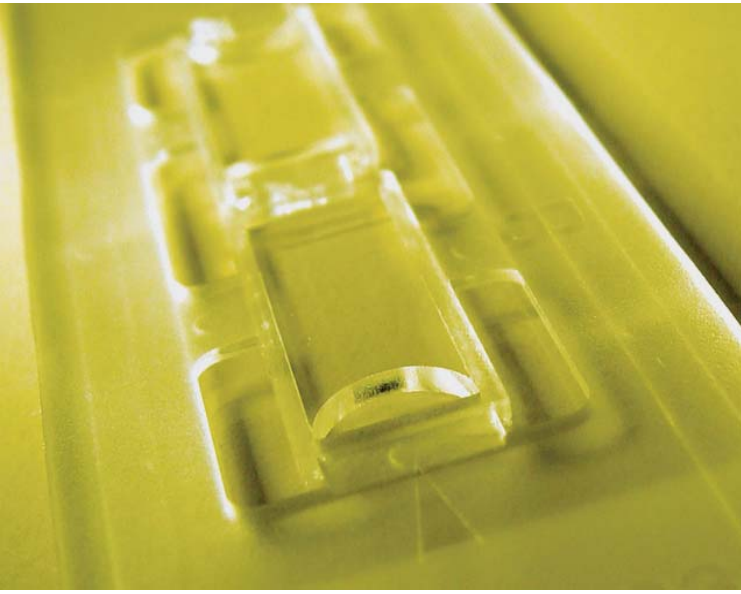

C-chip is a disposable plastic hemocytometer, solving the problems of conventional glass hemocytometer.

- Guarantees the repeatability and reliability by the complete quality control
- Two closed exact 100  $\mu\text{m}$  height chambers does not need a coverslip
- Eliminates the tedious washing step
- Reduces the exposure to infectious sample and the contamination
- Not fragile but very right and strong plastic

## Exact volume control & Easy to use

### Models

No.	Grid	Depth
DHC-N01	Neubauer Improved	100 $\mu\text{m}$
DHC-N02	Neubauer	100 $\mu\text{m}$
DHC-F01	Fuchs Rosenthal	200 $\mu\text{m}$
DHC-S01	Semen Test	10 $\mu\text{m}$
DHC-S02	No grid for CASA	10 $\mu\text{m}$
DHC-B01	Burker	100 $\mu\text{m}$
DHC-B02	Burker Turk	100 $\mu\text{m}$
DHC-T01	Thoma	100 $\mu\text{m}$
DHC-T02	Thoma New	100 $\mu\text{m}$
DHC-M01	Malassez	200 $\mu\text{m}$
DHC-P01	Petroff Hauser	10 $\mu\text{m}$

### Applications

Blood analysis (Hematology) :

- Blood cell counting

Cell Culture :



- Cell concentration measurement
- Cell viability test

Microbiology :

- Bacteria & Fungal spores counting

IVF, IUI :

- Sperm counting

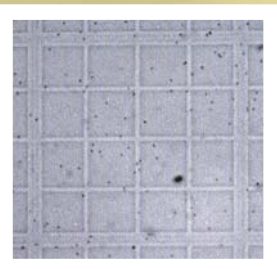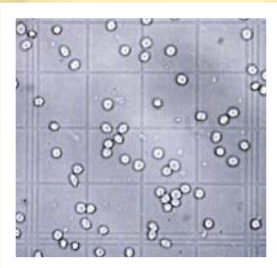

### Specifications

Dimension : 25 (W) x 75 (L) x Thickness 1.6 mm

Chamber depth : 100, 200  $\mu\text{m}$  (On sale)

10  $\mu\text{m}$  (Coming soon)

Chamber volume : 10 ~ 15  $\mu\text{L}$

### Ordering Information

#### Worldwide

Kyobo Tower, 1303-22, Seocho4-dong, Seocho-gu, Seoul, 137-070, Korea  
Tel : +82-2-3787-1964~5 Fax : +82-2-537-2649 E-mail : sales@incyto.com

#### Domestic (Korea)

Rm 511, SKC Institute, 911 Chongja-1 dong, Changan-gu, Suwon, 440-301, Korea  
Tel : +82-31-240-0331 Fax : +82-31-240-0312 E-mail : sales@incyto.com